
CORP CARE Miracle Smile

Gebrauchsanweisung.....	05
Instructions for Use.....	06
Mode d'emploi.....	08
Manuale di utilizzazione	09
Gebruiksaanwijzing	11
Használati utasítás.....	13
Návod k použití.....	14
Návod na použitie.....	16
Instrucțiuni de utilizare	17
Instrukcja obsługi.....	18

#A

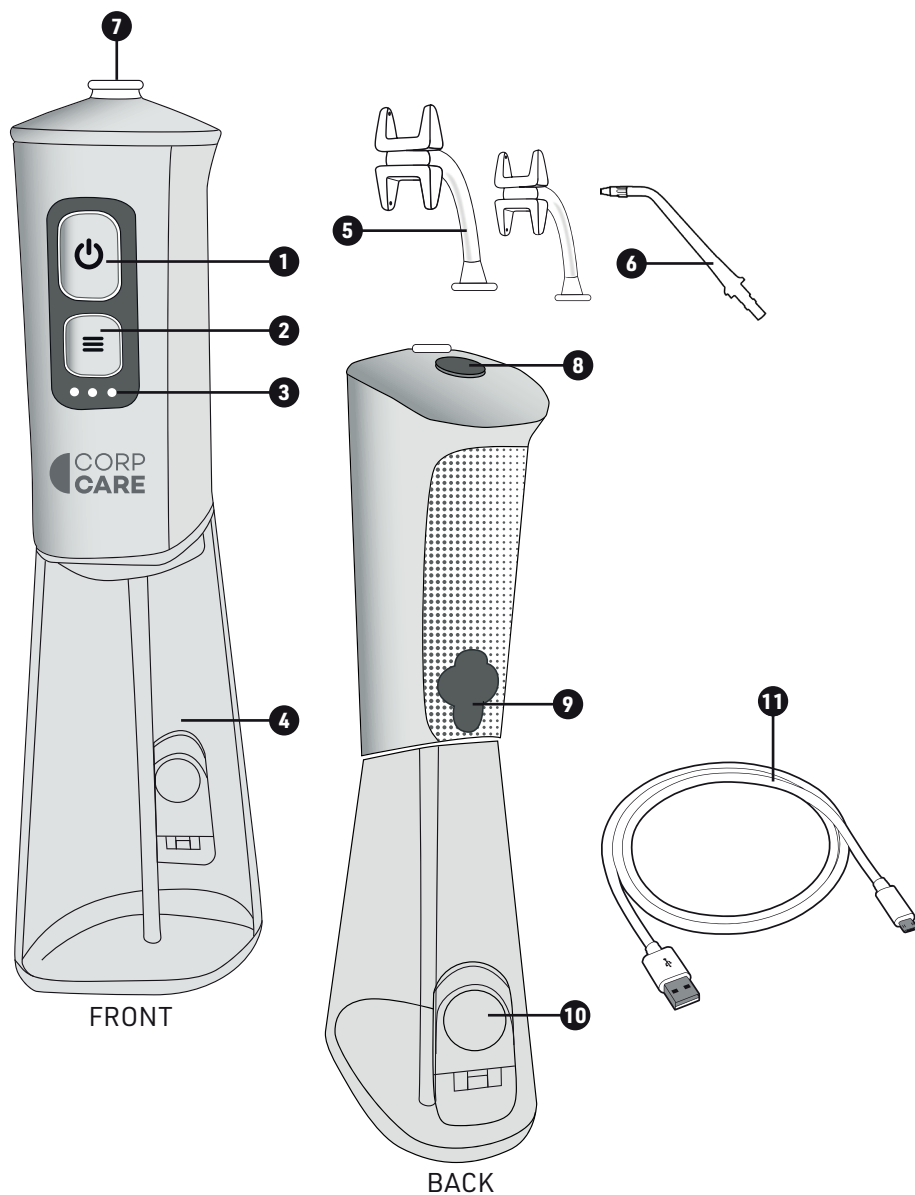DE

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Hergestellt in China

EN

IMPORTANT NOTES+SAFETY INSTRUCTIONS

- Please read the instruction manual carefully before using the appliance.
- This appliance is only suitable for oral hygiene.
- Keep the appliance out of the reach of children.
- This appliance may only be used by children and persons with reduced physical, sensory or mental capacities or lack of experience and/or knowledge if they have been supervised or instructed in the safe use of the appliance and can understand the hazards arising therefrom. Cleaning and user maintenance are not to be carried out by children unless they are supervised. Children must never be allowed to play with the appliance.
- Follow the instructions, use the appliance only as intended and/or as recommended by your dentist.
- If you experience discomfort or pain in your mouth, stop using the product and contact a dentist.
- Regularly clean the appliance according to the instructions.
- Check the appliance for damage or defects before each use, and do not use it if it is damaged.
- Do not use the appliance while it is connected to the power supply.
- If you have any questions or concerns, consult your dentist.
- Do not attempt to repair the appliance yourself! Contact customer service
- Do not direct the water jet directly under the tongue, or into the ear, nose or other sensitive

areas. This product can generate pressure that can cause serious damage to the areas to which it is applied.

- Accessories not included in the scope of delivery must not be used.
- Do not allow foreign objects to fall into any of the product's openings.
- Do not use iodine, saline or water-insoluble concentrated essential oils in this product. The use of such substances may affect its performance and shorten the product life.
- Only fill the container with warm water, mouthwash or other solutions recommended by your dentist.
- Remove any oral jewellery before using this product.
- Do not use the product if you have an open wound on your tongue or in your mouth.
- If you are planning a dental/oral hygiene procedure, consult your dentist before use.

CHARGING:

We recommend fully charging the appliance before first use.

- a. Open the protective silicone cover on the charging port on the back of the appliance **(#a)**.
- b. Plug the small end of the charging cable into the charging port and the large end into a USB connection or USB power adapter (not included) **(#b)**.
- c. The indicator lights will flash during charging. The appliance is fully charged when all three lights stop flashing and remain solid. Disconnect the charging cable from the appliance **(#c)**.

IMPORTANT - PLEASE READ BEFORE USE

Do not switch the appliance on until a cleaning attachment is attached to the appliance and this cleaning attachment is in your mouth! Always start by cleaning with the lowest intensity before slowly increasing it.

USE (#1 - #8)

CorpCare Miracle Smile is the perfect addition to your daily brushing routine.

1. Insert the H-shaped cleaning attachment or the jet nozzle into the attachment opening at the top of the appliance, and firmly press the desired attachment until it clicks into place.
2. Open the filling opening on the back of the water tank.
3. Carefully fill the water tank with water until it is full (200ml). TIP: For an enhanced experience, add an amount of mouthwash as desired before fully filling.
4. Fully close the filling opening. Refill the water tank as necessary.
5. Insert the attachment into your mouth. Place the lower part of the H-shaped cleaning head over the lower teeth - slowly close the mouth to position the upper teeth into the upper part of the H-shaped cleaning head.
6. After inserting it into your mouth, press the ON/OFF button. Use the mode button to select between low, medium and high pressure by pressing the button repeatedly; the indicator lights will light up accordingly.
7. Gently glide the head from side to side over your teeth in a continuous, controlled motion. If necessary, refill the water tank and repeat the process to achieve the desired result. **INFO:** A cleaning cycle takes about 10 seconds; after this, simply refill the tank. You can repeat the process as many times as you like. When you are finished, turn off the appliance by pressing the power button.
8. To remove the attachment, press the release button and carefully pull the attachment away from the appliance.

INFO:

The jet nozzle attachment is specially designed for the targeted cleaning of hard-to-reach areas, such as between teeth and below the gum line. Do not use the jet nozzle attachment on sensitive or inflamed areas. Hold the nozzle close to the tooth or desired position, move the appliance slightly to reach all sides of the tooth.

CLEANING + CARE OF THE APPLIANCE:

After each use, empty the water tank, remove any residue, and allow the appliance to air dry. Clean the

attachments thoroughly under warm water after each use. Store in a cool, dry place until next use.

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EN: This symbol means that such electrical equipment or batteries must not be disposed of with normal household waste. You are legally obliged to hand in such items - free of charge - to a public collection point or distributor-created collection point for the recycling of electrical devices and expired bulbs and batteries within the scope of the Electrical and Electronic Equipment Act that are not encapsulated by or fixed into the appliance, and can be removed or separated from the appliance non-destructively, for disposal at the end of their service life. Use rechargeable batteries instead of single-use batteries wherever possible. The recovery and recycling of old appliances makes an important contribution to protecting our environment. If disposed of improperly, toxic substances can enter the environment and have harmful effects on humans, animals and plants. Batteries with an increased pollutant content are also marked with the following symbols: Cd = cadmium, Hg = mercury, Pb = lead. Consumers are themselves responsible for deleting any personal data from appliances. Here are some best practices and recommendations to extend the life of your batteries and enable their reuse:

- Charge your batteries properly and fully to maximize their lifespan. If necessary, fully discharge the battery using a suitable appliance before recharging it.
- Consider the possibility of using rechargeable batteries. These can be reused and recharged multiple times, reducing the quantity of battery waste.

- Always use the correct types of batteries for your devices. Improper use may shorten battery life and potentially give rise to harmful effects.

Always follow the safety instructions when handling used batteries. There is a high risk of fire with old batteries containing lithium (Li). Therefore, special attention must be paid to the correct disposal of used batteries containing lithium. Incorrect disposal can also lead to internal and external short circuits due to thermal effects (heat) or mechanical damage. A short circuit can lead to a fire or explosion, and have serious consequences for both people and the environment. To avoid a short circuit, cover the poles or contacts of the battery before disposal.

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EU: MediaShop GmbH | Schneiderstraße 1, Top 1 | 2620 Neunkirchen | Austria
CH: MediaShop Schweiz AG | Leuholz 14 | 8855 Wangen | Switzerland