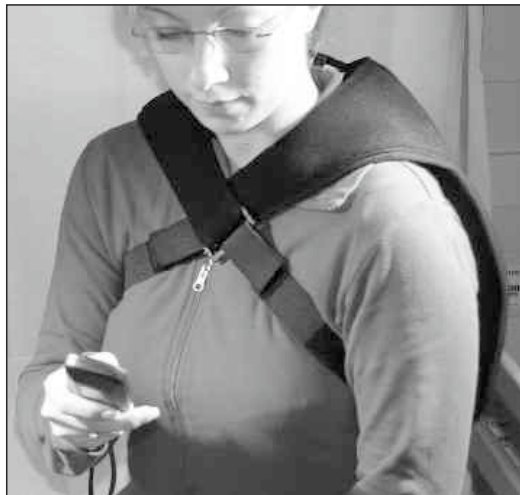
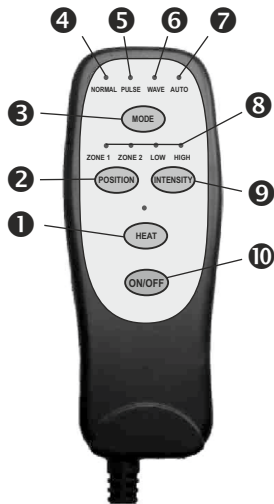
D Nacken- & Schultermasser

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GB Back & Shoulder-Massager

Important Safety Instructions

Read all instructions before using!

When using electrical appliances, especially when children are present, basic safety precautions should always be followed, including the following:

KEEP AWAY FROM WATER!

DANGER! - As with most electrical appliances, electrical parts are electrically live even when the switch is off. To avoid any risk:

1. Always "unplug it" immediately after use.
2. Do not use while bathing.
3. Do not place or store appliance where it can fall or be pulled into a tub or sink.
4. Do not use near, place in or drop into water or other liquid.

WARNING!

1. An electrical product should never be left unattended when plugged in and in operation.
2. Close supervision is necessary when this product is used by, on or near children or invalids.
3. Use this product only for its intended use as described in this leaflet. Do not use attachments not recommended by the manufacturer.
4. Never operate this product if it has a damaged cord or plug, if it is not working properly or if it has been dropped into water.
5. Keep the cord and product away from heated surfaces. Do not wrap the cord around the product.
6. Never use while sleepy or about to go to bed.
7. Never drop or insert any object into any opening.
8. Do not use outdoors or where moisture is evident.

SPECIAL INSTRUCTIONS

We recommend not to use an extension cord. If an extension is absolutely necessary,

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a 10-ampere cord should be used. Cords rated for less amperage may overheat. Care should be taken to arrange the cord so that it cannot be pulled or tripped over.

USING YOUR BACK & SHOULDER MASSAGER

Always have control pad in hand while use. Adjust the control to your comfort level.

1. Outfit yourself with the Back & Shoulder-Massager in the appropriate place as shown in the picture.
2. Connect the cord of the Massager to the AC Adapter, then plug into an outlet.
3. Select the desire massage mode/intensity and positioning combination.
4. You may use the Massager either standing, sitting or laying on a couch.

IMPORTANT WARNINGS

1. If you have any concerns about your health contact your practitioner, before using the device.
2. Heart pacemaker patients should not use this device.
3. In the unlikely event that you feel pain of the joints or muscles while using the device, it is recommended to interrupt the massage.
4. This massage unit is not to be considered a medical device. It is to simply massage your back and shoulder.
5. If you are currently undergoing medical treatment, it is strongly advised that you consult your physician before use.
6. Always start with lowest power setting.

The Back & Shoulder Massager has a built in thermostat which enables it to automatically shut off after 15 minutes.

REMOTE CONTROL FUNKTIONS:

HEAT: Press once to activate the heat function, or twice to switch it off.

POSITION: You can select a specific neck or shoulder region

ON/OFF SWITCH

INTENSITY: You may choose between two intensity levels: high, low

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PROGRAM SELECTION: Press the button to select one of four different massaging functions.

- ① Heat
- ② Position
- ③ Programselection
- ④ Normal
- ⑤ Pulsemassage
- ⑥ Wavemassage
- ⑦ Automatic program (run through all 3 programs)
- ⑧ Display of selected settings
- ⑨ Intensity
- ⑩ On/Off

Technical data:

Supply voltage: 12V

Current: 500mA

Input voltage: 100-240V

Adapter current: 0.2A

Dimensions: 38 x 47cm

Length of cable between adapter and device: 1.8 m

Operating temperature: -20°C - 60°C

Batteries, rechargeable batteries and electrical devices must not be disposed of with the normal household waste! All consumers are legally obliged to submit all batteries, rechargeable batteries and electrical devices, whether they contain pollutants or not, to a collection point in their community or district or to return them to a dealer, so that they may be disposed of an environmentally compatible way. Only submit batteries and rechargeable batteries in a discharged state!

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Maximex GmbH & Co. KG