Insektenstich-Helfer

Art.-Nr. 26102

- (DE) Gebrauchsanweisung**
- (GB) Instruction manual
- (FR) Notice d'utilisation
- (ES) Manual de instrucciones
- (IT) Istruzioni per l'uso
- (NL) Instructies
- (CN) 产品说明

INSEKTENSTICH-HELFER

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DE

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Art.-Nr. 26102

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Art.-No. 26102

Dear client,

congratulations on your purchase of the Insect bite pen. Please read these operating instructions thoroughly before use. It is important for proper use of the product.

1. How it works

When you are bitten or stung by an insect, protein-containing poisons are introduced into your body. The body then has an allergic reaction to these poisons. This involves the production of histamine, which causes itching (which can last for some time) as well as burning and swelling. To counteract this effect, the insect bite pen uses the principle of hyperthermia, or the controlled application of locally concentrated heat. When the temperature exceeds around 50°C, the protein chains of the poison are destroyed and the production of further histamine by the body is prevented. This principle is not new. Sailors have long used heated irons to relieve itching and the pain of bites and stings. The effectiveness of this approach has been clinically proven. Following immediate application, 74 % of test subjects experienced substantial relief from itching, burning and pain. 66 %


found that after 24 hours the sting had cleared up completely.

2. Technical data

- Gold-plated heating plates
- Temperature of heating plates controlled automatically at slightly over 50 °C
- Illuminated function display
- 9 V battery included
- Belt clip
- Pocket-sized
- Effectiveness clinically proven

Product description:

- a. Heating plate
- b. Control button with indicator lamp
- c. Battery compartment cover

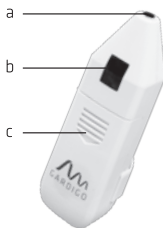3. Operation

Inserting Battery

To open the device, slide back the battery cover. A 9 V block battery is provided in the battery compartment on delivery. Take it out and unpack it. Connect the contact clip at the end of the cable to the battery, making sure that the battery polarity is correct. Once the battery has been connected, put it back into the battery compartment. The contact clip should point backwards, with the projecting cables stored underneath the device. When closing the battery cover, make sure that no cables are caught as they might thereby become damaged.

Changing the battery

If the indicator lamp does not light up any more when the control button is pressed, you have to replace the battery. Use a standard 9 V block battery for this. The battery provided will last for approx. 100 applications. Never use rechargeable batteries. Used batteries do not belong in domestic garbage. You are legally obliged to dispose of the batteries. Dispose of them through your specialist electrical supplier or at your local recycling point. Note: Batteries containing pollutants are marked with the following symbols: Pb = Battery contains lead, Cd = Battery contains cadmium, Hg = Battery contains mercury

4. How to use

Ensure that the insect sting is no longer in your skin. Place the heating plates of the bite pen on the area to be treated. Press the button and keep it pressed for the duration of the treatment time. The light will indicate that the heating plates are heating up. It takes about 25 seconds for the device to reach the required temperature of slightly over 50°C. Hold the device on the area affected for as long as you can comfortably bear.

Important: Children max. 25 seconds (no self-use of child)! Adults max. 30 seconds!

Attention! The heating mechanism does not switch-off itself automatically! If the insect bite is still itching after application, let the appliance cool down completely and repeat the application!

5. Information about the device

- This unit is intended for use on human beings.
- The device should be stored in such way as to prevent any inadvertent release of the button. Occasionally burns can occur as well as a severe build-up of heat.

- Only press the control button once the heating plate is in contact with the skin. If you press the button too early and “heat up” the heating plate, this can cause burns on the skin.
- For reasons of hygiene, the insect bite healer should only be used on one person.
- Diabetics should not use the insect bite healer as their pain sensitivity might possibly be reduced and they might inflict burns on themselves.
- This device should not be used with persons who are helpless, with children under 3 or persons insensitive to heat (e.g. persons with skin changes related to illness or persons with scarred skin in the area of application, after taking any pain-relieving medication, alcohol or drugs).
- The insect bite healer does not help against insect allergies.
- This unit is only intended for personal use and not for commercial applications.
- You should protect the device from impact, moisture, dust, chemicals, major temperature fluctuations and direct exposure to the sun’s rays.

6. Cleaning and disposal

- Always remove the battery from the device before you clean it.
- Clean the device with a slightly damp cloth. In case of more extreme soiling you can also moisten the cloth with mild soapy water.
- Make sure that no water penetrates into the device. If this should ever happen then only use the device once it is completely dry.
- Do not use any chemical cleaners or abrasive agents for cleaning.
- Please dispose of the device in accordance with the directive 2002/96/EG – WEEE (Waste Electrical and Electronic Equipment).

If you have any queries, please refer to the local authorities responsible for waste disposal.

Your Gardigo-Team

Environmental protection notice

At the end of its life span this product may not be disposed as normal household waste but must be disposed of at a collecting place for recycling of electrical and electronic equipment. The icon on the product, in the manual, or on the packing points to this fact. The materials are recyclable according to their labelling. Through reusing, material recycling, or other forms of utilisation of old devices you make an important contribution to the protection of the environment. Please ask the local administration for the responsible waste disposal centre.

Do not throw batteries or rechargeable batteries into household waste!

Consumers are legally obligated to return used and rechargeable batteries, whether they contain harmful substances* or not, to designated recycling areas, disposal sites, or stores where batteries/rechargeable batteries can be bought. By doing so, you are fulfilling all legal battery recycling and disposal obligations and are contributing to a better environment. Make sure batteries are completely discharged before disposing!

*marked with: Cd = cadmium, Hg = mercuric, Pb = lead

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