

ORAL IRRIGATOR



220/ 240V 50Hz 18W

INSTRUCTION MANUAL



- This appliance can be used by children aged from 8 years and above if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Cleaning and user maintenance shall not be made by children unless they are older than 8 and supervised.
- Keep the appliance and its cord out of reach of children aged less than 8 years.
- Children shall not play with the appliance.
- The appliance can be used by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and if they understand the hazards involved.

Please read the instruction manual before using the device for the first time and keep it for further reference.

## **IMPORTANT SAFETY INSTRUCTIONS**

1. Before connecting the appliance, check that the supply voltage indicated on the rating plate of the appliance matches that of your installation.
2. This appliance is only intended for indoor and domestic use.
3. Do not use this appliance for any other purpose than that described in this instruction manual.
4. Only put the appliance on a stable surface.
5. The appliance should not be left unattended when connected.
6. The appliance should not be installed in such a way that it might fall into the water or on the floor. The electrical outlet must be accessible at any time.
7. Never disconnect or handle the plug with wet hands.
8. Never move or disconnect the appliance by pulling on the power cord. Never put or store it in places where it might fall into a bathtub or shower.
9. Do not use it while bathing. Never use it while sleeping or drowsing.
10. Do not immerse in water or any other liquids.
11. Do not pick up an appliance that has fallen into a liquid without unplugging it first.
12. To avoid any risks of injury, do not direct the jet under the tongue, into the ear, nose, eyes or other sensitive areas
13. Check if the power cord is damaged before first use and before each use.
14. Only use the accessories supplied with the product.
15. To avoid any risks of injury, do not fill the water container with hot or boiling water. Do not put any object in the water container.
16. Fill only the water container with water or oral solutions.
17. Do not let the power cord hang or contact hot surfaces.
18. Do not use the appliance in places or rooms where aerosol and oxygen sprays have been used.
19. The appliance must be unplugged:
  - if there is an anomaly during operation,
  - if the power cord or plug is damaged,
  - before each cleaning or maintenance,
  - before mounting or dismounting,
  - if it is not used,
  - after use..
20. Never attempt to modify or repair the appliance yourself. Please contact an authorized after-sales service.

### **Health Warnings**

- Do not use if you wear a piercing. Remove it before using the appliance.
- Do not use if you have an open wound on your tongue, your lip (s) or in your mouth.
- Consult your dentist before using this appliance in the following cases:
  - If you are under a treatment or oral care,
  - If you have undergone important dental, oral or gingival care in the past 2 months,
  - If you have dental problems, or if you suspect having dental problems (fillings, crowns, etc.).
  - If your gums bleed excessively after using this appliance or if the bleeding persists.
  - If you have health problems,
  - If you wear a pacemaker or any other implanted device.






## PRODUCT DESCRIPTION



1. Tips
2. Water container
3. On / Off switch
4. Water Pressure Setting knob
5. Lid / Storage compartment for Tips
6. 360 ° Rotating Tip
7. Tip eject button
8. Pause Button

Volume of Water Container: 600 ml

## Accessories

				
Ordinary High-Pressure Tip(x3)	Periodontal Tip (x1)	Orthodontic Tip (x1)	Dental Plaque Tip (x1)	Tongue Cleaner Tip (x1)

\* Ordinary High-Pressure tips: 3 colors of rings for easy identification of the user in a family for hygiene reasons.

## OPERATION

### Step 1

#### Filling the Water Container

1. Remove the lid and the water container from the base
2. Fill the water container with warm water. Then carefully wipe the water container.
3. Put the water container back on the base, with the graduation in front of you, and press firmly on it until you hear a click. Then put back the lid on the water container.




### Step 2

#### Inserting and removing tips

Insert the tip into the center of the hole at top of the dental jet handle. A click is heard, indicating that it has been correctly installed.

**Note:** To remove the tip from the dental jet handle, press the tip eject button and pull out the tip.

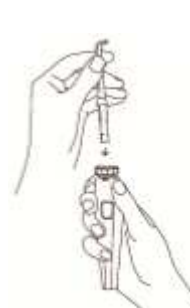

Inserting



Removing

### Step 3

#### Power ON

1. Plug the power cord into an electrical outlet. If the outlet is controlled by a wall switch, make sure it is turned on.
2. Lean over the sink and put the tip between your teeth and your gums.
3. Change the jet orientation by turning the knob at the top of the handle (see illustration below).
4. Make sure the water pressure setting knob is positioned before position 2. With the other hand, press the ON / OFF switch (position I).



### Step 4

#### Setting pressure parameters

1. Turn the water pressure setting knob to position 2 for first use.
2. Gradually increase the pressure until you reach the unit pressure that suits you.

**Note:** The jet pressure modes are shown below:

**Gentle Mode** : positions 1-2-3

**Intermediate Mode** : positions 4-5-6-7





**Strong Mode** : positions 8-9-10



### Step 5

#### Recommended cleaning technique

1. Adjust the jet at a 90-degree angle on your gums, close your lips a little to allow water to fall down into the sink and to avoid splashes. For best results, start with the area of the molars (back teeth) to gradually move towards the front teeth.
2. Continue to do this until all areas around and between teeth are cleaned.
3. You can temporarily stop the waterflow at any time by pressing and holding the pause button on the handle during the pause. Release the button to continue.
4. When finished, turn off the appliance by pressing the ON / OFF switch (position 0). Unplug the appliance.
5. Empty the water container completely to avoid the proliferation of bacteria in the stagnant water.

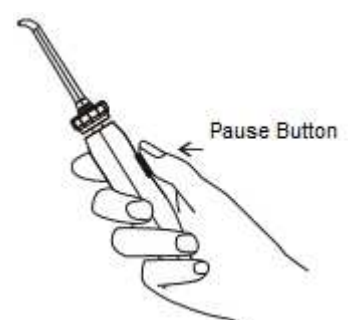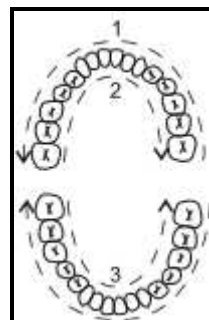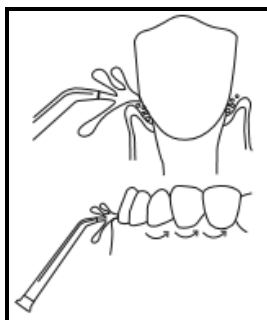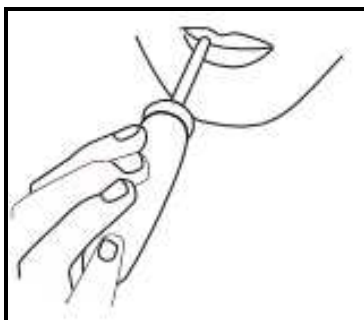


## PROCEDURES FOR A CORRECT USE OF THE TIPS

Before using the tips, you should lean over your sink and operate as described below according to the tip used.

### Periodontal tip

The periodontal tip propels a jet of water or oral solution deep into the periodontal pockets.

1. Put the soft tip against a tooth, under the gums at a 45-degree angle.
2. Turn the knob of pressure to position 2
3. Gently clean the gums and insert the tip between the teeth. Continue to follow along the gums until you finish cleaning.

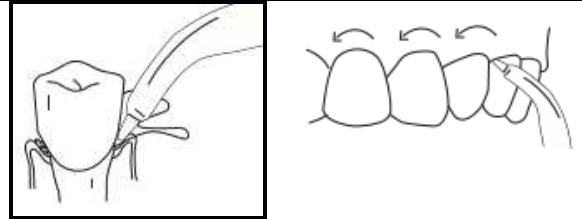

### Tongue Cleaner


1. Put the tongue cleaner on your tongue and turn on the appliance.
2. Select the pressure unit 2 and increase until you reach the unit of pressure that suits you.
3. Gently put the tongue cleaner in the center of your tongue and gently push forward.



### Orthodontic tip

The orthodontic tip is only designed to brush and rinse hard-to-reach areas simultaneously.

Use the brush to clean around braces or dental prosthesis. Water will flow through the brush to dislodge food residues.

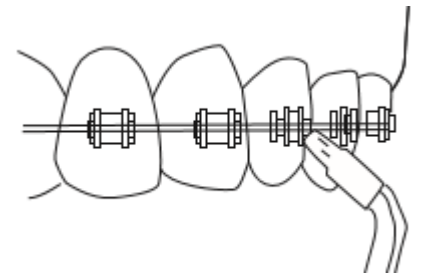

### Dental plaque tip

Direct the jet at a 90-degree angle towards your gums, close your lips a little to allow the water to fall down into the sink and to avoid splashes. Start from the molars (back teeth), slide gently along the gum and near the teeth (the bristles of the toothbrush should touch the teeth), and stop briefly between the teeth. Direct the spray gradually towards the front teeth.

Continue until all areas around and between the teeth, inside or outside, up and down are cleaned.



## CLEANING AND STORAGE



- Before cleaning the appliance, put the ON / OFF switch to the "O" position, rotate the setting knob before position 2 and disconnect the appliance from the electrical outlet.
- Do not immerse the appliance and the power cord in water or any other liquid.

### WATER CONTAINER

- After use, empty the water container and rinse it carefully if you have used an oral solution.
- To dry the water container, put it upside down on a flat and stable surface.
- To avoid ingestion, do not put isopropyl alcohol or any other cleaning liquids in the water container.
- In order to prolong the life of the oral irrigator and prevent it from becoming dirty, partially fill the water container with tap warm water from time to time. Direct the tip of the dental jet towards your sink. Turn on the appliance. Let it operate until the water container is empty.

### BODY OF APPLIANCE

- Do not use steel wool or scouring powders.
- if necessary, clean the body of the appliance with a soft, slightly wet cloth. Dry it thoroughly.

### **TIPS**

- Before removing the tip, empty the water contained in it by directing the tip towards the sink;
- Press the eject button of the tip. Wipe the handle and the outside of the tip. Let it dry.
  
- Keep the appliance indoors, in a dry and safe place.
- Keep the appliance away from sunlight, moisture and rain.
- Keep the appliance out of reach of children.

## **TROUBLESHOOTING**

<b>Problems</b>	<b>Causes</b>	<b>Solutions</b>
Dental jet tips and tongue cleaner: Water leaks from the tip connection and the dental jet handle.	Tip not completely connected.	Remove the tip and reinstall it.
Periodontal tip: damaged tip	The tip is worn.	Replace with a new tip.
Inadequate pressure	Water container not correctly positioned.	Adjust the position of the water container again. Press firmly on it.
The appliance does not work.	Engine is not running.	Please make sure that the plug works.

**WARNING: You should not dispose of this device with your household waste. A selective collection system for this type of product is implemented by your local authorities. Please contact your local authorities to find out how and where collection takes place. These restrictions apply because electrical and electronic devices contain dangerous substances that have harmful effects on the environment or on human health and must be recycled. This symbol indicates that electrical and electronic devices are collected selectively. The symbol shows a waste container crossed out with an X symbol.**

**If the power cable is damaged, it must be replaced by the manufacturer, his after sales service, or a person qualified in this area in order to avoid any accidents arising.**

Imported by PRODIS SAS, 1 rue de Rome, 93110 Rosny-sous-Bois, France/ Made in the P.R.C  
Marketed by EUROtops Versand GmbH, Elisabeth-Selbert-Str. 3 D-40764, LANGENFELD, GERMANY

