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**Warning:** Please read this instruction manual carefully before using the instrument and keep it for future use.

#### A. **SAFETY INSTRUCTIONS**

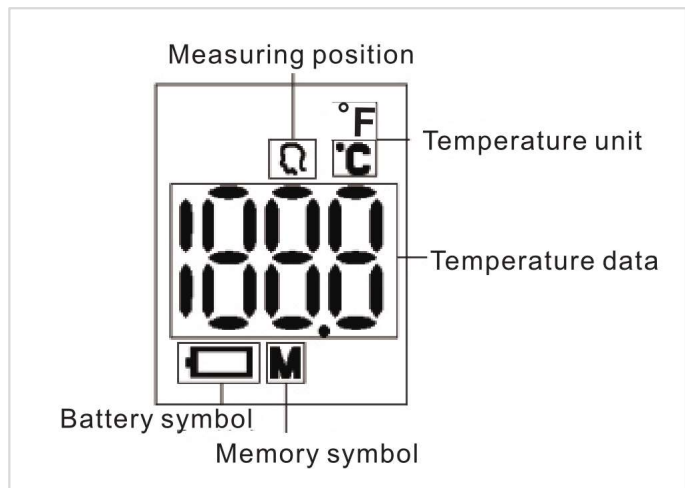
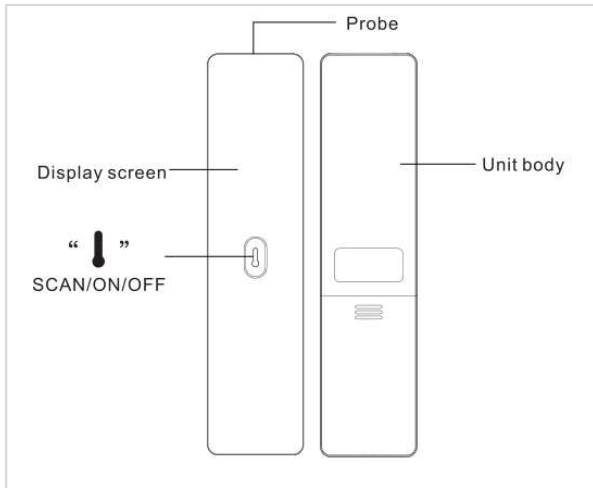
1. The device is intended exclusively for private and not for commercial use. Only use this device as describe in this instruction manual.
2. Do not place any part of this appliance in water or any other liquid.
3. Do not handle the unit or plug with wet hands.
4. Keep it away from heat sources such as radiators, to avoid deformation of plastic parts.
5. Do not use coarse or abrasive sponges/cloths or steel wool to clean the appliance.
6. This unit is not designed to be operated by means of an external timer or a separate remote-control system.
7. It is suggested to try the device by yourself first, so as to know how to use it.
8. In order to avoid affecting accuracy of the device due to the external too cold or too hot environment, it must be placed in general environment of room temperature of 16°C -35°C for at least 30 minutes before use to avoid incorrect measurement results.
9. The person measured should be kept at room temperature for at least 20 minutes to make body temperature balanced.
10. After strenuous exercise, you should rest for at least 30 minutes before taking measurement.
11. Wipe the skin dry and push aside hair strands before commencing.
12. When repeated measurements are made continuously, move away the probe between each measurement and rest for at least 5s before making the next to obtain the most correct result.
13. The normal reading value is about 35.5 °C ~37.8 °C . The user should take more temperature to know his or her own health temperature. Pay attention to physical changes or consult your doctor if there is any abnormality.
14. Check the probe is clean before each measurement to ensure accuracy.
15. Do not use the thermometer under temperature extremes (below 16°C/60.8°F or over 35°C/95°F) or humidity extremes (below 15%RH or over 80%RH).
16. Do not expose the thermometer to temperature extremes (below -20°C/-4°F or over 55°C/131°F) or humidity extremes (below 15%RH or over 93%RH).
17. Don't touch or blow infrared sensor. A polluted or broken infrared sensor may cause inaccuracy.
18. Do not use a mobile phone nearby when the thermometer is working. Portable and mobile RF communications equipment can affect MEDICAL ELECTRICAL EQUIPMENT. Falling to do so may cause inaccuracy.
19. Do not modify this equipment without authorization of the manufacturer. The device is not repairable and contains no user serviceable parts.
20. The operator shall not touch the thermometer and the patient simultaneously.

#### B. **PRECAUTIONS D'USAGE**

- Always place the thermometer in the same place on your forehead when taking a forehead temperature, as the value may be different from one place to another.
- In order to ensure the accuracy of the measurement, there should be no hair, sweat, cosmetics and dirt, etc. on the forehead when measuring.
- If you place ice or the like on the forehead of the feverish patient the measurement result will be falsified. Users should avoid action in this case.

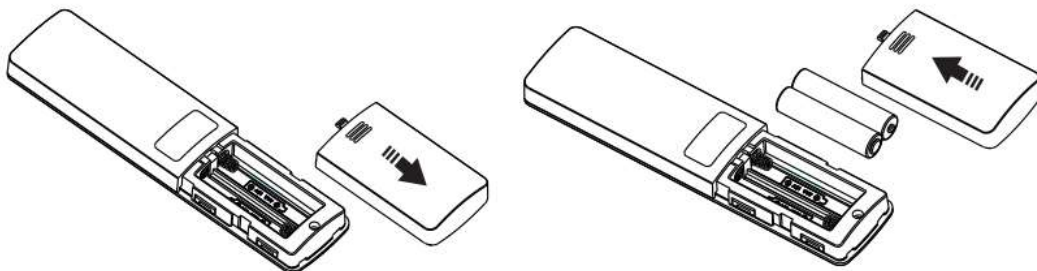
**C. INTRODUCTION**

**D. HOW TO USE IT**

**1) Insert or replace batteries**

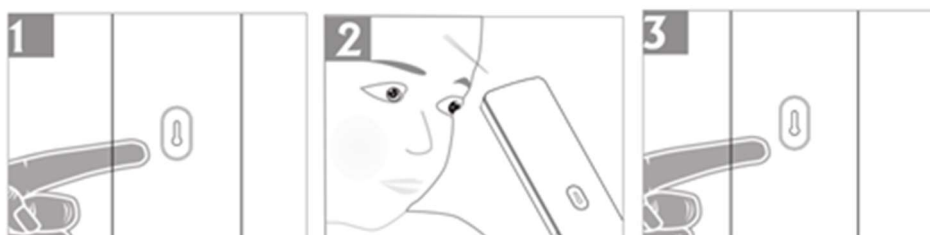
- Remove the battery cover on the back of the unit.
- Insert 2 AAA batteries observing the correct polarity indicated.
- Replace the cover.
- When the battery icon “” is indication, please replace the old battery immediately
- Please take out the battery if you do not use the thermometer for a long time.

**2) Instruction for use**

Average temperature of human body is around 35,5 °C ~ 37,8 °C. Consult a doctor if needed.

- Turn the thermometer on by pressing “” switch.
- A beep sound will be then heard and the screen will shows °C ---- . The unit is now ready for use.
- Place the thermometer’s probe point to the forehead (distance 15mm to 50mm), press “” and release after ~ 1 second, the thermometer will beep and display the temperature. If a consecutive reading is required then . repeat the above please wait For at least 5 second procedure as often as required.

- Location indication: The icon indication the correct measuring position when measuring.
- Temperature data indication: The temperature reading is displayed after measuring.

- Press the “” switch for 5 seconds to turn off the thermometer. Or the unit will switch off automatically after 60 seconds.

- **If a reading below 32 °C is recorded a “Lo” sign will be displayed followed by 2 consecutive beeps.**
- **If a reading above 37.8°C is recorded there will be 6 consecutive beeps.**
- **If a reading above 42.2°C is recorded a “Hi” sign will be displayed followed by 2 consecutive beeps.**

### 3) Recalling Readings from Memory

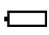
- Ensure the thermometer is switched off.
- Press the “” button for 4 seconds to put the device on memory mode. The last reading will be displayed.
- Press and release the “” button again to display the next reading in a sequence from the last reading to the first reading.
- The thermometer has the memory to recall the last 32 readings.
- The thermometer will automatically turn off if it idles for more than 20 seconds.

Sometimes parents tend to rely on their principal in demand on their child's forehead, it is even a reflex whether the child is feverish or not. Or, it suffices that this hand is cold or hot to misjudge the temperature. With the forehead thermometer, the results obtained are quite precise, provided you follow the correct procedure.

**You will add 0.2 ° C to the evaluated value to get the exact internal temperature.**

It is dangerous for patients to perform a self-evaluation and self-treatment based on the measuring results. Be sure to follow doctors' instruction. A self-evaluation may cause deterioration of diseases.

## E. TROUBLESHOOTING

| Error message  | Problem   | Solution  |
|--|---|---|
| Hi   | Temperature taken is higher than 42.2°C(108°F).           | Operate the thermometer only between the specified temperature ranges. In the event of a repeated error message, contact your retailer or Customer Services |
| Lo   | Temperature taken is lower than 32°C(89.6°F).             | Operate the thermometer only between the specified temperature ranges. In the event of a repeated error message, contact your retailer or Customer Services |
| Err  | Sensor not welded   | Contact your retailer or Customer Services  |
| ErE  | EEPROM is abnormal  | Contact your retailer or Customer Services  |
|  | Low power   | Replace the batteries   |
| ErH  | The temperature of the measuring environment is too high. | Decrease the temperature of the environment, keep it within 16°C-35°C.  |
| ErL  | The temperature of the measuring environment is too low.  | Increase the temperature of the environment, keep it within 16°C-35°C.  |

## F. CARE AND MAINTENANCE

### 1) Maintenance

## T 5001 / AET-R1F1

Remove any stains on the body with a soft and dry cloth.

Gently wipe the dirt on the body with a soft cloth, use a cotton swab or soft cloth stained with alcohol, and gently wipe the lens of the temperature measuring head. After using the product, the product can be sterilized with 75% medical alcohol on the soft towel or cotton sliver.

Don't wash with water or detergent containing abrasive or benzene.

Do not immerse in liquids.

### 2) Storage

The device must not be stored or used at an excessively high or low temperature or humidity, in sunlight, in association with an electrical current or in dusty locations. Otherwise may occur inaccuracies.

Remove the battery if non-use for long time.

### 3) Calibration

This thermometer is calibrated at the time of manufactured. If the thermometer is used according to the instruction, periodic recalibration is not required. If at any time you question the accuracy of the measurement, please contact the retailer immediately. Don't attempt to modify or reassemble the thermometer.

## G. SPECIFICATION

- Material: ABS
- Power: 2 x AAA battery
- Resolution 0.1
- Operating Conditions: 5°C-40°C (41°F to 104°F) with a relative humidity of 15%-95%
- Storage Conditions: -20°C-+55°C (-4°F to 131°F) with a relative humidity of 15%-95%
- Forehead temperature: Measurement range: 32.0°C-42.2°C
- Measurement accuracy:  $\pm 0.2^{\circ}\text{C}$

*The European directive 2012/19/EC about electronic and electric waste, requires that you can't throw away defective domestic appliance with common waste. Used device has to be collect separately to optimize the recuperation rate and the recycling of the materials to reduce the impact on the health and the environment.*

*This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capacities, or by persons lacking in experience or knowledge, unless they are supervised or instructed on the use of this apparatus by a person responsible for their safety. Children should be supervised to ensure they do not play with this appliance.*

### **Precision and precautions about the battery use**

- *Throwing away a used battery in the environment or with household waste pollutes and avoids the recuperation of reusable materials.*
- *It is important to limit the battery consumption and to respect the following rules:*
  - *Favor alkaline battery (last longer than saline battery) and when it's possible use rechargeable battery.*
  - *Leave used batteries and accumulators in adapted containers that you can find in some retailers. In that way, metals don't pollute the environment and the health.*
- *Battery has to be inserted by respecting the polarity show on the apparel and the battery. An incorrect positioning can damage the apparel, cause leaking and even provoke fire and the explosion of the battery.*

T 5001 / AET-R1F1

- *To assure a good functioning, battery has to be in a good condition. If there is a problem, put new battery.*
- *Never try to recharge battery not made for this. It could leak, warm up, cause fire or explode.*
- *Replace all the battery at the same time. Never mix up alkaline battery with saline battery.*
- *Used battery has to be removed.*
- *Also, remove the battery if you're not using the apparel for a long time.*
- *Never try to short-circuit the terminal of a battery.*
- *Never throw away the battery in fire because they can explode.*
- *Take out the accumulator from the apparel before reload it.*
- *We recommend that an adult watch if a child change the battery to respect those rules.*
- *If you swallow a battery, go immediately to your doctor or the hospital.*
- *For the best ongoing performance, we recommend using only fresh, high-quality, alkaline batteries.*

**Importé par / Importiert von / Geïmporteerd door / Imported by :**

**CAPTELEC, 59170 Croix - FRANCE**

Carton et emballage papier à  
trier / Kartons und  
papierverpackung recyceln /  
Verpakkingsmateriaal in karton  
en papier te sorteren / Recycle  
carton and paper packaging