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HAPPY FEET STIMULATOR MASSAGER

EN

Warning : Please read the instruction carefully before using the product and using for future use.

NOTES BEFORE USE

1. People with the following conditions should consult their physician before using the product:
 - Patients who are weakened by a severe disease
 - Who uses a medical electronic device, which is embedded in the body, such as a pacemaker.
 - Who suffers from heart disease, hypertension
 - Who suffers from backbone injury or disease.
 - Who is suffering from bone fragility, especially on the backbone.
 - Who has some abnormality or deformity of limbs or backbone.
 - After surgery.
 - Who is pregnant or in a period of monthly illness.
 - Who suffers from malignant tumors, cancer, hyperthyroidism, tuberculosis, advanced diabetes
 - Who is suffering from osteoporosis.
 - Who is suffering from blood circulation problems or thrombosis.
 - Who is suffering from skin disease or infection.
 - Who has skin issues or a fragile skin
 - Who has skin allergies
 - Who is undergoing medical treatment or feels some physical abnormality previously or is being treated now.
 - Who has injuries, diseases due to illness or accident.
2. Product is not intended for use by the physically or mentally impaired (including children).
3. Children should not use this product. They should be under supervision to ensure that they will not treat this product as a toy.
4. Do not apply medical or cosmetic creams on the product or on your foot before using this product.

SAFETY PRECAUTIONS

1. Only use the product in accordance with the instructions in this user manual.
2. Use the product only with the accessories provided. Do not connect any other accessory.
3. Recommended usage time for an individual is 25 minutes at a stretch.
4. Make sure that your skin is clear and clean before placing the electro pads.
5. Don't place the electro pads on the chest

S 3003 / AST-300H

6. Do not place the electrodes near the chest, which could increase the risk of heart fibrillation
7. Do not let the electro pads in contact with metallic substance like bracelet, watch, necklace etc... when the device is on.
8. Do not let water enter the product or be splashed near the surface of the product.
9. Stop using the device immediately if you feel nausea or uncomfortable when using the appliance.
10. Don't sit or stand on the device. Don't throw it away and don't let it fall.
11. Water contact with electricity is dangerous. Do not leave the product in a humid environment (such as in the bathroom or shower room or swimming pool).
12. Before connecting the product to the power supply, check that the product voltage and supply voltage are the same.
13. Do not use the product if cover is damaged or if there is a mechanical or electronic failure.
14. In order to prevent the risk of electric shock, do not open the product outer cover.
15. If the product or the power cord is damaged, do not try to repair the product. This must be handled by a specialist / technician otherwise the product guarantee will become void.
16. If an unusual sound or noise is being heard from the product, turn off the product and immediately disconnect the power connection. Please contact an authorized OTO Service center immediately.
17. Do not use the product if it is damaged or if it has a mechanical or electronic failure.
18. Do not use this product in high temperature conditions or beside hot objects such as heaters, stoves.
19. Do not destroy, bend stretch, twist or cut the power cord to avoid electric shock.
20. Make sure to use this device in a distance of 1 meter minimum away from other devices transmitting microwaves or electromagnetic in high frequency or short wave.

PRODUCT DESCRIPTION

- 1 Foot massager
- 2 Adaptor 100-240 V
- 3 2 output wires
- 4 4 electrode pads
- 5 Remote control

CONTROL PANEL

A. Modes

Modes	EFFECTS
1	Acupuncture pushing
2	Acupuncture
3	Acupuncture kneading
4	Acupuncture taping
5	Scrapping
6	Squeezing
7	Massage
8	Pushing massage
9	Pushing squeezing
10	Acupuncture squeezing
11	Acupuncture hammering
12	Kneading
13	Thumping
14	Scrapping pressing
15	Cupping
16	Body shaping
17	Hammering
18	Massage tapping
19	Pushing
20	Rolling pounding
21	Squeezing
22	Stroke
23	Acupuncture therapy massage
24	Shiatsu
25	Rolling kneading


B. Sole

99 intensity level


C. Body (Electrode pads)

99 intensity levels

CONTROL BUTTONS

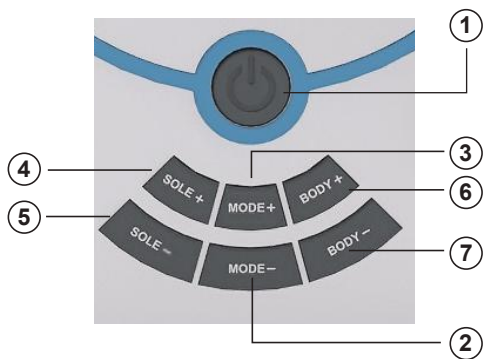- 1 ON/OFF
- 2 Choose one mode - downward
- 3 Choose one mode - upward
- 4 Increase the output intensity of foot sole
- 5 Decrease the output intensity of foot sole
- 6 Increase the output intensity of electrode pads
- 7 Decrease the output intensity of electrode pads

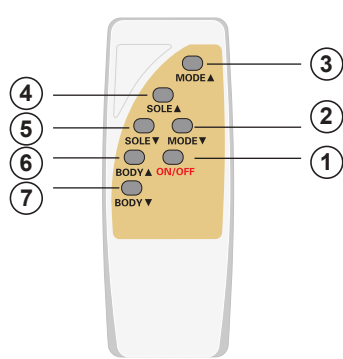CONNEXION OF THE ADAPTOR

Connect the output wire of the adaptor in the whole on the left of the foot massager, then plug the adaptor. When plugging the adaptor, a blue light will come out the LCD screen very shortly.

CONNEXION OF ELECTROPADS

Connect each output wire of the electrode pads [3] (please refer to the product description).
 Take off the protection layer of the electro pads
 Connect the output wires of the electro pads on the right of the device.
 Place the two electro pads on the body parts you want: calves, thighs, waist, butt, lower back, back

Key functions:

To operate the remote control, put two AAA batteries (not provided) in the battery compartment in correct polarities.

Buttons	Functions
MODE	Please check correctly that the device is well plugged. Press the button ON/OFF during 2 seconds on the device or on the remote control, a blue light will come out the LCD screen showing the different settings. Place your feet on the sole place. Place the electrode pads on the body part you want (except for the chest) Select the stimulation mode and intensities that you desire.
STOP	Press the button ON/OFF anytime. If there is no weight on the device or if the electrodes are not placed on the body, the device will turn off automatically.
INTENSITE	99 intensity levels If you don't feel anything, do not hesitate to gradually increase the intensity levels.
CHALEUR	No heat comes out the device.

Auto timer: 25 minutes. Power will turn off automatically after 25 minutes.
 Always unplug the AC power cord after use.

To really feel the effect, place the two electrode pads on your own body. If the two electrode pads are placed on two different bodies, the massage from the electrode pads won't work.

You can roll your foot sole on the mechanical roll with round picks [6] (please refer to the product description) for a little extra to your stimulation massage.

Every time when changing the modes, you will have to reprogram the intensity.

TROUBLE-SHOOTING

This page summarizes the most common problems you could encounter with the product :

Problem	Possible Cause	Suggestion
The product does not work	The product is not plugged in properly.	Ensure the voltage on the product corresponds to the local mains voltage and then plug in again.
	The electrode pads are sliding or fell during the session	Remove them, clean them with a damp rag and replace them.
	The product is likely to be overused.	Unplug the product and let it cool down before using it again.
	The intensity is weak	Increase slightly the intensities to be more convenient for you.
	A heavy charge was on the product	Reduce the surplus of the charge and turn it on again.
The skin gets reddish and irritated.	The electro pads are stick tightly on the skin.	Remove them and place them gently on your skin again.
	The electro pads are dirty	Clean them with a damp rag
The device turns off during the massage session	The device is not well plugged	Unplugg, plug it again and turn it on.
	The massage session is over	After 25 minutes, the device will turn off automatically

CARE AND MAINTENANCE**Cleaning**

Do not use steel wool, abrasive cleaning agents or corrosive liquids (such as gasoline or acetone) to clean the product.

1. Unplug the plug from the power outlet.
2. Wipe with a dry cloth to clean the product.
3. Do not use a damp cloth or detergent

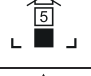
Storage

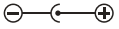
1. Unplug the plug from the power outlet.
2. Store the product in a cool dry place away from sunlight.

3. Do not put a heavy weight load on the product.

If the cover is damaged or there is any mechanical failure, please stop using the product.

Explanation of figures, symbols, warning statements and abbreviations on the equipment

No.	Symbol	Description
1		BATCH CODE. This symbol should be accompanied by batch code.
2		Symbol for "SERIAL NUMBER". This symbol shall be accompanied by the manufacturer's serial number.
3		DATE OF MANUFACTURE. This symbol shall be accompanied by the date of manufacture.
4		CONSULT INSTRUCTIONS FOR USE.
5		Symbol for "MANUFACTURER". This symbol shall be accompanied by the name and the address of the manufacturer.
6		KEEP DRY.
7		Type BF Applied Part.
8		FRAGILE.
9		THIS SIDE UP.
10		DO NOT STACKING.
11		Symbol for "ATTENTION, CONSULT ACCOMPANYING DOCUMENTS".

12		Class II Equipment.
13		Indoor use only.
14		RoHS Compliant Product.
15		Polarity of d.c. power connector.
16		The product and the accessories are compliant to the european directive

Specifications:

Mode NO.: AST-300H

Rated: 5 V dc, 1000 mA

Powered by adaptor: EE-0501000DPE

Rated Voltage: 100-240V Rated frequency: 50-60Hz

Rated input: 0.1 A Rated output: 5V 1000mA

Output waveform: Square wave Pulse duration: 115 µs

Pulse repetition frequencies: 10-100HZ Effect of load impedance: 1KΩ

Maximum amplitude of output voltage: 72V

Software version No:MC0188F-REV-V1.0

Operation condition:

10 – 40 °C

30% RH ~ 75%RH

860 hPa to 1060 hPa

Store and transport condition:

0 – 40 °C

≤80% RH

860 hPa to 1060 hPa

The European directive 2012/19/EC about electronic and electric waste, requires that you can't throw away defective domestic appliance with common waste. Used device has to be collect

separately to optimize the recuperation rate and the recycling of the materials to reduce the impact on the health and the environment.

This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capacities, or by persons lacking in experience or knowledge, unless they are supervised or instructed on the use of this apparatus by a person responsible for their safety. Children should be supervised to ensure they do not play with this appliance.

Importé par / Importiert von / Geïmporteerd door / Imported by :
CAPTELEC, 59170 Croix - FRANCE

Model number: AST-300H

This equipment belongs to Type BF of Class II according to Clause 5.

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Fabricant/Hersteller/Fabrikant/Manufacturer

SHENZHEN OSTO TECHNOLOGY COMPANY LIMITED

ADD: No.43, Longfeng Road, Xinsheng Community,

Longgang Street, Longgang District, 518116 Shenzhen

City, Guangdong Province, China.

Version:A.2

2018-1-19

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Share Info GmbH