POWER GROW COMB

by Velform®

 Instruction Manual

 Mode d'emploi

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 Manual de instrucciones

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ISL GmbH
Karl-Schiller-Strasse 4a
26209 Hatten-Sandkrug Germany
www.isl-de.com

IMPORTANT INSTRUCTIONS RETAIN FOR FUTURE USE.

SAFETY INSTRUCTIONS

Read all the instructions carefully for safe use of the Power Grow Comb.

Never remove the warning label from the electrical device.
KEEP THESE INSTRUCTIONS IN A SAFE PLACE.

* Power Grow Comb by Velform is distributed with a warning label in English. (to know where to find this label, see "Inserting the batteries" in the instructions section of this manual).

* Before using Power Grow Comb for the first time, stick the attached warning label in your language on top of the warning label in English.

* Do not direct the laser beam towards people or animals. Do not look directly into the laser beam even from a great distance. This device generates class 2M laser radiation in accordance with European standard EN 60825-1. It is possible to inadvertently blind other people with this type of laser.

* To ensure safe operation of the device, always have the device repaired by a professional using original parts. Do not allow children to use the device unattended, as they could inadvertently blind other people with it.

* Never use this device as a substitute for medical care.

* If you suffer or have suffered from a serious illness, consult your doctor before using Power Grow Comb by Velform.

* This device should not be used for purposes other than those described in the manual.

* Discontinue use if pain or swelling occurs after application.

* Discontinue use immediately if any kind of irritation occurs.

* Do not use on any part of the body not indicated in the instructions for use.

* This device is not recommended for use by people under 18 or by people with disabilities without adult supervision.

* This device should not be used by pregnant women.

* Do not use this product while asleep.

* Power Grow Comb by Velform is not a medical or therapeutic product.

* Do not dismantle the product under any circumstances. It should only be dismantled or repaired by authorized professionals.
* This product should not be used by people (including children) with physical, sensory or mental disabilities, or if the lack experience or knowledge about it, except if they have been supervised or instructed about the safe usage of the product by a responsible person. Children should always be supervised to make sure they don't play with the product.

The power Grow comb by Velform® includes –

A - The power Grow Comb: Hair growth unit.

B - Massage comb: Relaxing rubber massage comb to use in the shower.

C - Manicure set: Luxury pampering travel kit.

D - Hair brush: Professional hair salon brush.

E - Instructions

Introduction

Congratulations you have purchased the most advanced hair growth comb on the market today.

Please read and follow the instructions carefully.

Information about the Power Grow Comb by Velform®

The Power Grow Comb by Velform® uses Lasers and infrared LED lights to help stimulate scalp cells. This ground breaking therapy also helps improve blood circulation. The combination of Lasers and infrared lights helps to nourish hair follicles.


The Power Grow Comb by Velform® also stimulates hair follicles by vibration technology. The patented vibration therapy will help increase blood flow to your hair follicles and in turn help improve hair growth.

The attached comb also enables the Laser and light therapy to reach your hair follicles. The teeth of the combs are positioned before or after the Laser and LED's to create perfect hair separation.

Through scientific investigation, engineering and design the Power Grow Comb by Velform® will help grow your hair!
INSTRUCTIONS


Inserting the Batteries

1. Open the battery cover.
2. Insert 4 AAA size batteries as shown:

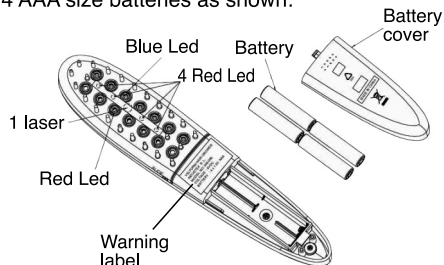Turning the Power Grow Comb by Velform® on.

To turn on the unit press the power switch by the logo. You will hear a short beep and the LED by the power switch will glow red. This means the Power Grow Comb by Velform® is on.

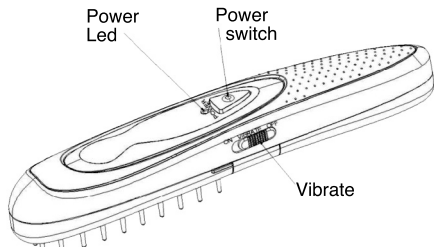
Turning on the vibration mode.

The Power Grow Comb by Velform® also features a vibrate mode. On the side of unit, move the VIBRATE switch to the 'ON' position. The vibration mode will then be activated. Note – the vibration mode will only work when the unit is on and the power button has been pressed.

Use Power Grow Comb according to your own schedule and convenience: during 5 or 10 minutes a day, or during 10 to 15 minutes every other day. You will begin to notice results in approximately 8-12 weeks.

Detaching the Hair Comb

If you desire, it is possible to detach the comb of the Power Grow Comb by Velform® and apply the laser and infrared lights directly over your scalp. To detach the comb please locate the 'slide' labels at either side of the comb. Once located please slide the comb away from the unit. To re-attach the comb please click the unit back into its original position. You will hear a 'clicking' sound.

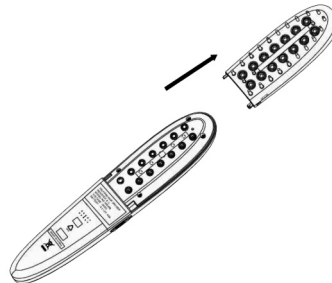Directions of Use

You should use the Power Grow Comb by Velform® without applying any lotions or gels.

Use Power Grow Comb according to your own schedule and convenience: during 5 or 10 minutes a day, or during 10 to 15 minutes every other day. You will begin to notice results in approximately 8-12 weeks.

Please move the unit slowly around your head in the direction of hair growth, from front to back, and bottom to top. Please make sure you cover all of your head and move in a slow controlled motion.

Massage Comb


The Massage Comb is especially designed to help circulate the blood flow in your scalp and improve your hair growth.

It is possible to use the Massage Comb with either wet or dry hair.

Comb in the shower to help apply shampoo.

Please place your fingers in the holes of the Massage Comb. (See Figure)

Move the massage comb around your scalp in a circular and controlled motion.

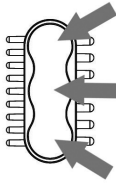SPECIFICATION

Power Grow Comb by Velform® (model n°: WB/096)

Input voltage

6V DV (4 x AAA size battery)

Laser wavelength

650nm +/- 5nm

Laser class

0,93 mW Class 2M certified

Red LED wavelength

660 nm +/- 5nm

Blue LED wavelength

470 nm +/- 5nm

Output

420mW without vibration

600m W with vibration

Power Grow Comb by Velform® complies with the following directives:

* Safety of Laser Products: IEC 60825-1 / A2 IEC 60825-1

* Electromagnetic Compatibility: EN 61000-6-1 (2007) + EN 61000-6-3 (2007)

* Restriction of Hazardous Substances (RoHS)

* The product has been granted the CE mark for its compliance with the above directives.

WARRANTY

This product is covered by a warranty against manufacturing defects subject to the time periods stipulated in the legislation in force in each country.

ATTENTION

DISPOSING OF MATERIALS

The symbol of a crossed out rubbish bin with wheels indicates that you must be aware of and follow local regulations with regard to the disposal of this kind of product.

Do not dispose of this product the way you would dispose of normal household waste.

You must follow the relevant local regulations. Electric and electronic products contain hazardous substances that have harmful effects on the environment and human health and should be recycled properly.

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